

HYDRAULIC EMERGENCY STRETCHER



■ FEATURES

EMERGENCY STRETCHER

■ BACK-RAISE

Can be adjusted with the handle gas spring.

■ SIDE GUARD

Can be used for transfer the patient to bed.

■ MATTRESS BASE

Compatible mattress bases are used for X-Ray.

■ BUMPER

Available at every corner to reduce impact deflection.

■ I.V POLE STORAGE

■ KNEE-RAISE

Can be set with using the crank.

■ CENTRAL LOCK CASTOR

Equipped with an 8" Castor with a central locking system, there are pedals in all four corners for easy access.

■ OXYGEN TUBE HOLDER

Located under the mattress, with a capacity of 1m³.

■ HYDRAULIC PEDALS

Set the height position trendelenburg and anti-trendelenburg.



■ FEATURES

EMERGENCY STRETCHER



5TH WHEEL

Equipped with a 5th wheel for product mobility more stable.



X- RAY CASSETTE HOLDER

Provided for use of X-Ray Casette for the imaging process by providing a clear internal structure of the patient's body.



ANGLE INDICATOR

Assist medical staff in monitoring the patient's position, and ensuring the patient is in a safe and comfortable condition.



HANDLE KNEE-RAISE

Helps to organize knee-raise.



POST HANDLE

Provided with a handle that can be opened and closed, to help push the stretcher.



DOWN ADJUSTMENT PEDAL

Helps to lower the base of the mattress properly simultaneously. It's in the foot position.



CASTOR

Equipped with 4 castors measuring 8 inches with a central lock system.



I.V POLE

There are 4 hooks and can be placed in 4 positions.


EMERGENCY STRETCHER

31214




- Emergency Stretcher with side guard features and hydraulic mechanism, and equipped with 5th wheels to make it easier to maneuver and control the stretcher when used.


DIMENSION (mm) ±

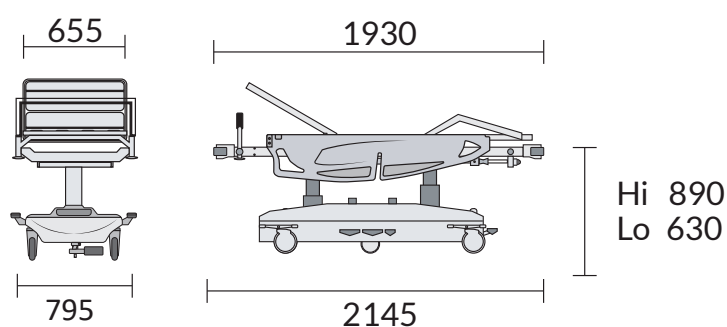



STRETCHER POSTITIONS

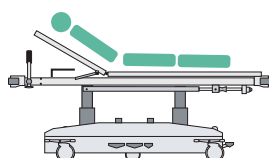

Back-raise angle can be adjusted with a slope of 0° - 70°.

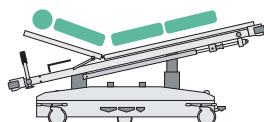

Trendelenburg - Anti - Trendelenburg angle can be adjusted with a tilt of ± 15° - 15°.

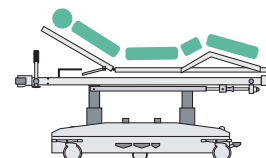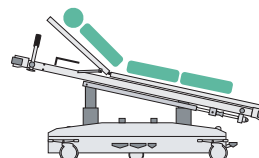

Knee-raise can be adjusted up to ± 0 - 25°.

■ TECHNICAL SPECIFICATIONS

EMERGENCY STRETCHER

Construction	Made from steel pipes and steel sheets processed with the latest manufacturing technology.
Mattress Base	Made from materials that can be used for X-Ray procedures, with a frame quality steel.
Side Guards	Made from polypropylene material. It can be adjusted in position and is equipped with a lock.
Position	<ul style="list-style-type: none"> • Mattress base height (low-high) \pm 630 - 890 mm - hydraulic mechanism • Trendelenburg tilt angle - anti trendelenburg \pm -15° - 15° - hydraulic mechanism • Backrest tilt angle \pm 0 -70° - gas spring mechanism. • Footrest tilt angle \pm 0 -25° crank mechanism.
Mattress	Made of foam covered with synthetic leather.
Castor	4 single wheel 8" castors with central lock system (one step on all four castors locked). 1 castor as a fifth wheel, its use can be adjusted to facilitate mobility.
Accessories	Stainless steel infusion pole and belt strap.
Finishing	Epoxy powder coating (oven paint).
Max Load	170 kg.
Safe Workload	200 kg.
Product Weight (\pm)	135 kg.
Dimensions (\pm)	\pm (H) 2145 x (W) 795 x (H) 920 mm - 1180 mm (overall). \pm (H) 1930 x (W) 655 x (H) 630 mm - 890 mm (at mattress base).

Specifications are subject to change without prior notice





EMERGENCY STRETCHER

31214



EMERGENCY STRETCHER with Folding Table

31215



Certified Company
Cert No.: 16 00 L 15191



Certified Company
Cert No.: 01 05 L 15006



Certified Company
Cert No.: 08 104 0021



Certified Company
Cert No.: 03 141 20026

PT MEGA ANDALAN KALASAN

MARKETING OFFICE

Gran Rubina Business Park Lt. 22C,
Komp. Rasuna Epicentrum,
Jl. HR. Rasuna Said, Jakarta,
Indonesia 12940

t. (+62) 21 - 837 00 555
f. (+62) 21 - 837 00 335
marketing@mak-techno.com

International Contact:
export@mak-techno.com

FACTORY

Jl. Tanjung Tirto 34,
Tirtomartani, Kalasan, Km 13, Yogyakarta
Indonesia 55571

t. (+62) 274 - 497 069
f. (+62) 274 - 496 226

www.mak-techno.com